

# ENDOMAX<sup>®</sup> PLUS

## Cordless Endodontic Handpiece



- Cordless convenience offers complete portability.
- Easy to read LCD display.
- Heavy-duty high torque motor.
- Automatic forward/reverse functions.
- Switch from rotary to reciprocating with a press of a button.
- Five torque settings and nine rotation speeds available for different applications.
- High capacity lithium-ion battery.

# ENDOMAX<sup>®</sup> PLUS

## Cordless Endodontic Handpiece



Activation Button

Mode Selection Button


Torque Limit Button

RPM Setting Button

Up/Down Selection Buttons


Power On/Off Button

### Auto-Reverse Forward:

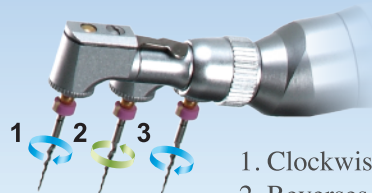


1. Clockwise rotation.
2. Reverses when torque limit is reached.
3. When the load is removed, rotation continues forward.

### Reciprocating Mode:

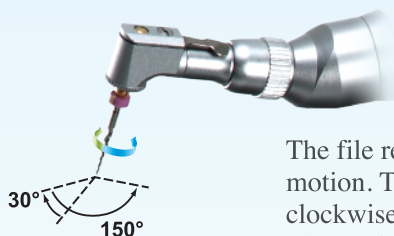

The file reciprocates in a back and forth motion. The default setting is 30° clockwise and 150° counterclockwise. The angles can be adjusted in 20° increments from 30° to 170°.

### Auto-Reverse Stop:

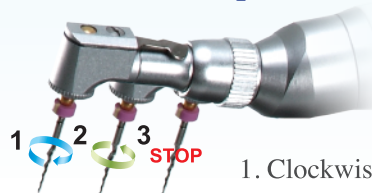

1. Clockwise rotation.
2. Reverses when torque limit is reached.
3. When the load is removed, rotation stops.

### SPECIFICATIONS

INPUT	AC 100V~240V, 50/60Hz
FUSE	1A 250V
BATTERY TYPE	3.6V 2040mAh
DIMENSIONS (HANDPIECE)	L201 xW26 xH29 mm
DIMENSIONS (CHARGING STATION )	L105 xW85 xH59 mm
WEIGHT (HANDPIECE)	155g
WEIGHT (CHARGING STATION )	125g

## DENTAMERICA<sup>®</sup>

18688 E. San Jose Ave.  
City of Industry, CA 91748 U.S.A.  
Phone: (626) 912-1388  
Fax: (626) 913-0510  
www.dentamerica.com  
www.dentamerica.asia

DISTRIBUTOR